


2021□□□□□□□□□□□□□□□□4□□□□□□□□□□

ICARE2 - FM10

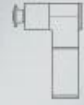

# FASCIAL MASSAGER

Boots your energy everyday



Gear 1  
Muscle wake



Gear 2  
Fascia relaxation



Gear 3 lactic  
acid decomposition



Gear 4  
Professional mode



# Small size and powerful


MAGIC GENIUS

Muscle stretching/massage relaxation



4 Gears

Relax muscles and  
relieve pain

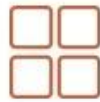

4 Pro massage head

Precise massage, relax  
whole body



Smart Timer set

Science guard  
human health




≈0.5KG

The weight of a bottle  
of mineral water



2500mAh

Equipped with large capacity  
aluminum battery

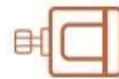

Low noise

Brushless motor  
stable & strong


\* Data from the plant laboratory, a slight error, please understand.

## Applicable scene

MAGIC GENIUS

---

It is suitable for a variety of people and occasions

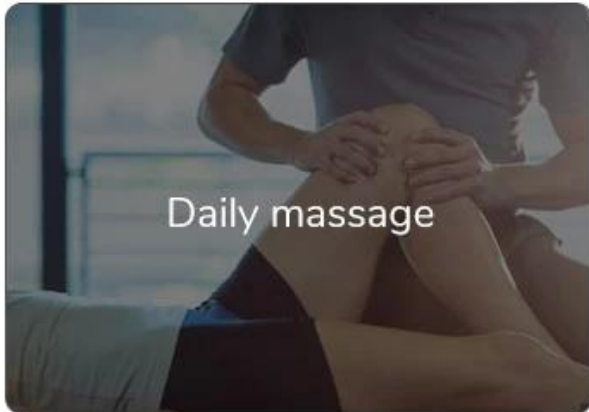

# How to use

Simplify to easy one-button control

MAGIC GENIUS

Switching indicator light

Long press 3s power on  
Long press 3s power off  
Single click to switch gear

• ON/OFF

• Charging socket

• Power light



## Four gear strength

MAGIC GENIUS

Daily body massage



Gear **1**: 1650rpm/min

Muscle wake



Gear **2**: 1800rpm/min

Fascia relaxation



Gear **3**: 2400rpm/min

lactic acid decomposition



Gear **4**: 3000rpm/min

Professional mode




The first time to use the suggestion from the lowest level,  
non-professionals do not use high gear!

# Professional massage head

MAGIC GENIUS


Matching soft massage head for different parts



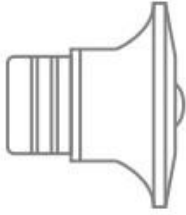


### Flat head

Suitable for each part of the muscle relaxation, shaping

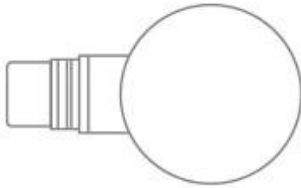


### Spherical head

Suitable for massaging arms, waist, back, buttocks, thighs, Large muscle groups such as the calves

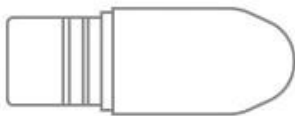

### Cylindrical head

Suitable for impacting deep tissues such as nerves, palms, soles, etc.

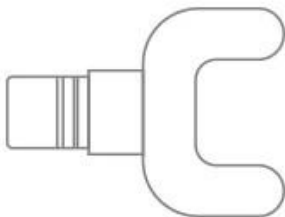

### U-shaped head

Suitable for massaging the neck, spine and Achilles tendon



# Professional massage

## Use as directed

MAGIC GENIUS

---



**Shoulder muscles**

5-20sx3times




**Hip muscles**

15-20sx3times




**Waist muscles**

10-60sx3times

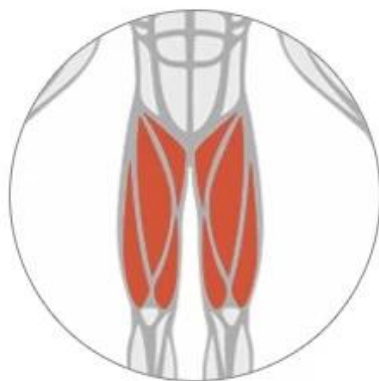

**Leg muscles**

15-45sx2-3times

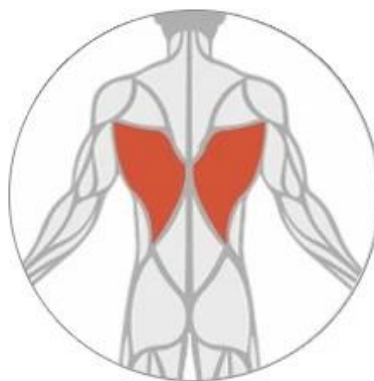

**Back muscles**

30-60sx3times

---

If you feel unwell, please stop immediately. For more usage methods, please read the instructions

---

# Product Parameter

MAGIC GENIUS



Brand: Magic Genius

Item No.: iCare2- FM10

Interface: USB Type C

Input Power: DC 5V

Battery capacity: 2500mAh

Rated power: 28.5W

Website: [www.Magic-Genius.com](http://www.Magic-Genius.com)